



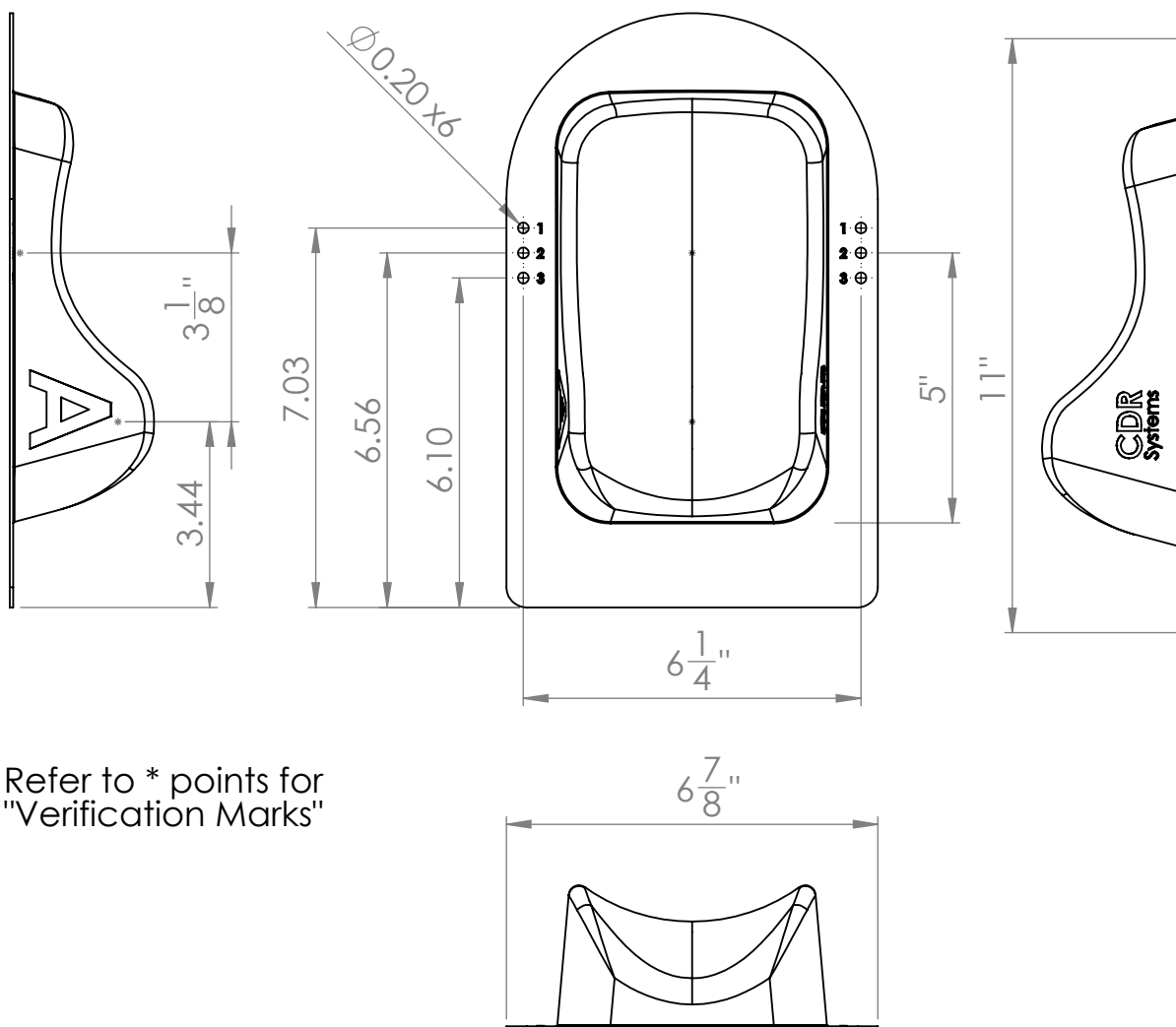

APPLICATION: RIGID MEDICAL POSITIONING HEADREST ACCESSORY FOR RADIOTHERAPY, MRI, AND CT.
ENVIRONMENT: 100% MR-SAFE AND RADIOLUCENT. METALLIC/MAGNETIC MATERIALS STRICTLY PROHIBITED.
STERILITY: SUPPLIED NON-STERILE. REUSABLE CLINICAL DEVICE.

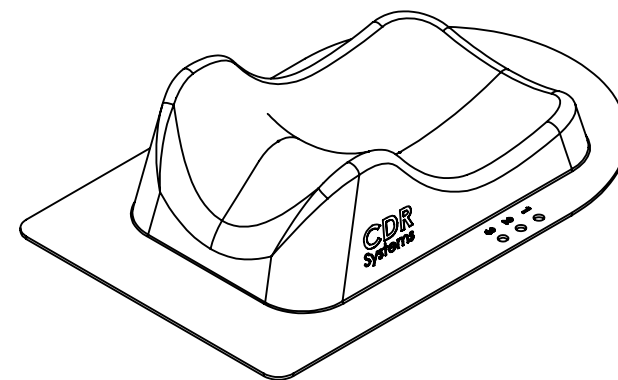

BASE MATERIAL: PETG SHEET.
 NOMINAL RAW SHEET THICKNESS: 1.5mm (0.060in)
BIOCOMPATIBILITY: MATERIAL MUST WITHSTAND REPEATED WIPEDOWNS WITH 70% ISOPROPYL ALCOHOL (IPA), QUATERNARY AMMONIUM DISINFECTANTS, AND BLEACH-BASED SANITIZERS WITHOUT CRAZING, CLOUDING, OR DEGRADATION.
OPTICAL/SURFACE: CLEAR WITH SMOOTH FINISH.
PRIMARY PROCESS: VACUUM THERMOFORMING OVER MALE TOOLING.
SECONDARY PROCESS: CNC ROUTING / TRIMMING FOR OUTER CONTOUR AND INDEXING HOLES.
SURFACE/EDGE FINISH: NO CHIPS, BURRS, WEBBING, DRAW LINES, OR TOOL MARK INCLUSIONS PERMITTED. DEBURR AND SMOOTH TO PREVENT SHARP EDGES OR PATIENT DISCOMFORT.
QUALITY SYSTEM: MANUFACTURING, PROCESSING, AND PACKAGING MUST COMPLY WITH AN ISO 13485 CERTIFIED QUALITY MANAGEMENT SYSTEM.

ADDITIONAL NOTES:

- SALES CODE: CHS26-A
- SEE 21-001-01 FILE FOR MANUFACTURING

UNLESS OTHERWISE NOTED ALL DIMENSIONS ARE INCHES TOLERANCES(Inch): ANGULAR: .X = ± .030" ± 0.5° .XX = ± .010" .XXX = ± .005" TOLERANCES(mm): X = ± 1.0 mm X.X = ± 0.25 mm X.XX = ± 0.13 mm PROPRIETARY AND CONFIDENTIAL THE INFORMATION CONTAINED WITHIN IS THE PROPERTY OF CDR SYSTEMS AND SHALL NOT BE REPRODUCED, DISTRIBUTED OR DISCLOSED TO OTHERS EXCEPT WHEN SPECIFICALLY AUTHORIZED IN WRITING BY CDR SYSTEMS.	DRAWN BY JUSTIN Y.			 CDRSystems Patient Positioning Systems with you in mind.
	APPROVED BY Jade Denis			
	APPROVED DATE (M/D/Y) 07/22/2026			
	SHEET 1/1			
	SIZE A			
REV 0			MANUFACTURING PROCESS VACUUM FORMING	DRAWING NUMBER 21-001-01

0	INITIAL RELEASE	07/22/2026	JY
REV	DESCRIPTION	DATE(M/D/Y)	NAME