



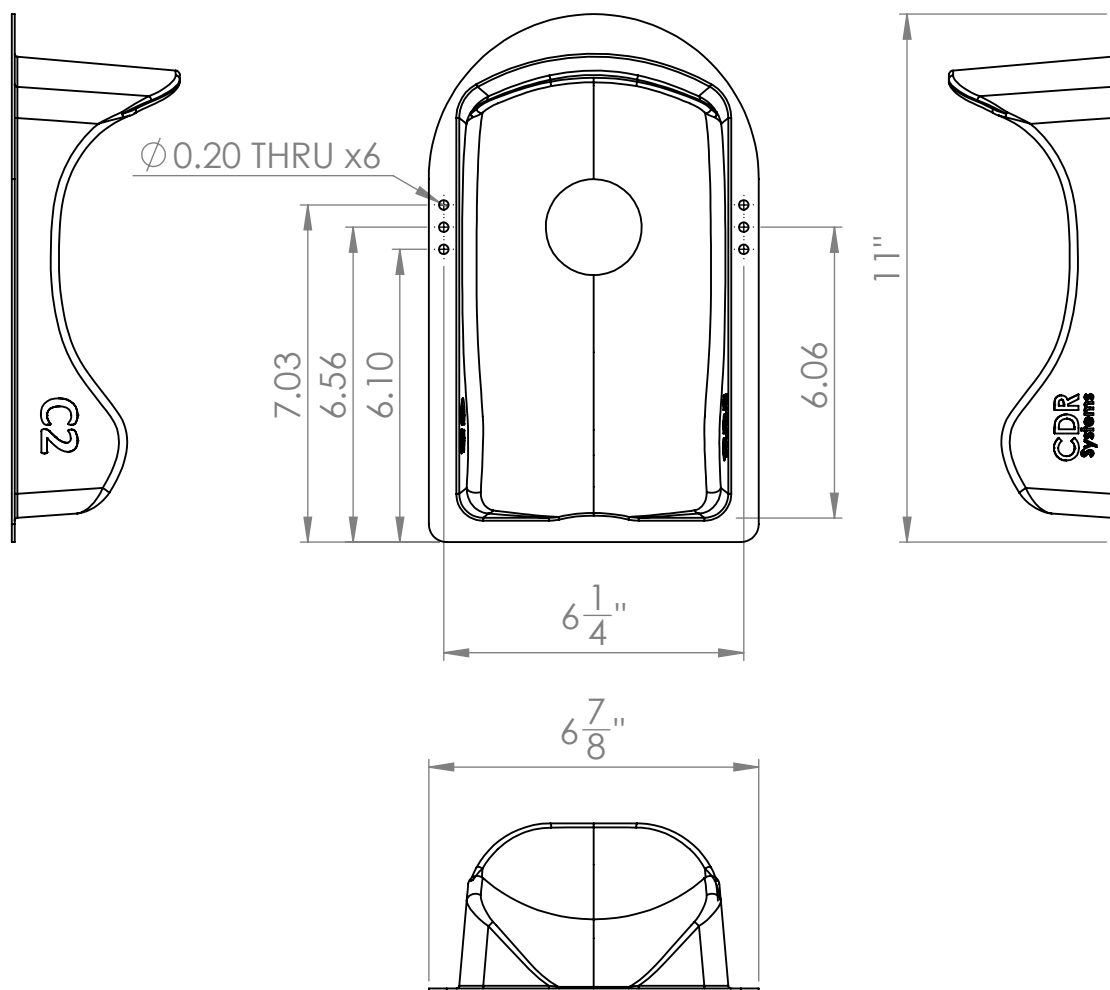

APPLICATION: RIGID MEDICAL POSITIONING HEADREST ACCESSORY FOR RADIOTHERAPY, MRI, AND CT.
ENVIRONMENT: 100% MR-SAFE AND RADIOLUCENT. METALLIC/MAGNETIC MATERIALS STRICTLY PROHIBITED.
STERILITY: SUPPLIED NON-STERILE. REUSABLE CLINICAL DEVICE.

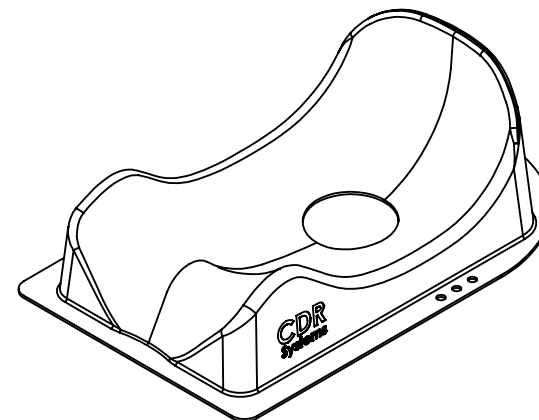

BASE MATERIAL: PETG SHEET.
 NOMINAL RAW SHEET THICKNESS: 1.5mm (0.060in)
BIOCOMPATIBILITY: MATERIAL MUST WITHSTAND REPEATED WIPEDOWNS WITH 70% ISOPROPYL ALCOHOL (IPA), QUATERNARY AMMONIUM DISINFECTANTS, AND BLEACH-BASED SANITIZERS WITHOUT CRAZING, CLOUDING, OR DEGRADATION.
OPTICAL/SURFACE: CLEAR WITH SMOOTH FINISH.
PRIMARY PROCESS: VACUUM THERMOFORMING OVER MALE TOOLING.
SECONDARY PROCESS: CNC ROUTING / TRIMMING FOR OUTER CONTOUR AND INDEXING HOLES.
SURFACE/EDGE FINISH: NO CHIPS, BURRS, WEBBING, DRAW LINES, OR TOOL MARK INCLUSIONS PERMITTED. DEBURR AND SMOOTH TO PREVENT SHARP EDGES OR PATIENT DISCOMFORT.
QUALITY SYSTEM: MANUFACTURING, PROCESSING, AND PACKAGING MUST COMPLY WITH AN ISO 13485 CERTIFIED QUALITY MANAGEMENT SYSTEM.

NOTES:

- SALES CODE: CHS26-C2
- SEE 21-004-02 FILE FOR MANUFACTURING

UNLESS OTHERWISE NOTED
 ALL DIMENSIONS ARE INCHES
 TOLERANCES(Inch): ANGULAR:
 .X = ± .030" ± 0.5°
 .XX = ± .010"
 .XXX = ± .005"
 TOLERANCES(mm):
 X = ± 1.0 mm
 X.X = ± 0.25 mm
 X.XX = ± 0.13 mm
 SURFACE FINISH 125
 PROPRIETARY AND CONFIDENTIAL
 THE INFORMATION CONTAINED WITHIN IS THE PROPERTY OF CDR SYSTEMS AND SHALL NOT BE REPRODUCED, DISTRIBUTED OR DISCLOSED TO OTHERS EXCEPT WHEN SPECIFICALLY AUTHORIZED IN WRITING BY CDR SYSTEMS.

DRAWN BY			 CDR Systems <i>» Patient Positioning Systems with you in mind.</i>		
JUSTIN Y.					
APPROVED BY					
-					
APPROVED DATE (M/D/Y)			C2 EXTENDED HEADREST		
08/06/2026					
SHEET	SIZE	REV	MANUFACTURING PROCESS		DRAWING NUMBER
	1/1	A	0	VACUUM FORMING	21-004-02

0	INITIAL RELEASE	08/06/2026	JY
REV	DESCRIPTION	DATE(M/D/Y)	NAME